CHAMBERS AND GLOVE BOXES

Touch screen, automatically

CONTROLLED ATMOSPHERE ISOLATOR/GLOVE BOX

The NEW **#890-THC** accurately maintains temperature and humidity (Rh) inside a roomy single operator glove box. We also offer multiple operator units.

The unit is operated by merely touching the color control screen for the functions desired. It includes your private password protected administrator controls. Touch screen viewing of the data logging function is included as well as serial port connection for data logging.

Automatic purging cycles can be programmed. The humidity range control accuracy is \pm 2% Rh.

Temperature can be regulated in either "F" or "C" degrees and the Peltier style thermoelectric assembly is accurate to \pm 2 degrees C. Temperature range is from 16 to 44 degrees C.

Features included are:

- Password protected, configurable "PID" controlled.
- Inert gas hookup capability.
- White ambidextrous Hypalon™ gloves.
- Platinum RTD with sealing gland.
- Four (4) pre-set conditions.
- 24 hour one touch, on screen data logging.
- Room temperature and Rh read outs.
- Thermoelectric assembly (Peltier device).
- Nebulizer for the Rh system.
- The top section is one piece clear acrylic plastic for easy viewing. It features "Easy Clean" corners.
- Four (4) key cock valves for purging. Two (2) on the main chamber and two (2) on the transfer chamber.
- 8" (203 mm) Color panel touch screen. Date and time display.
- Humidity (Rh) transmitter with 8 feet of sensor cable.
- Standard 24 month warranty (Not gloves or consumables).
- Serial port for external computer data logging & remote access.
- Transfer chamber is 11" (280 mm) I.D. x 12" (305 mm) long.
- "Bright Light" system (40,000 hour lamp guarantee).

Optional Drying Train package. Catalog #890-THC/DT

The Drying Train accessory includes two (2) canisters filled with Molecular SieveTM and circulation (vacuum) pump. This accessory can be used if a house supply of inert gas is not available. Periodic regeneration of the Molecular SieveTM is recommended for optimal drying performance. Refer to Accessories Section.

NOTE: If drying train feature is preferred, order as: #890-THC/DT (Domestic) or #890-THC/DT/EXP (Export).

"NEW PRODUCT"




#890-THC & #890-THC/EXP

Product Specifications

Inside Dimensions	41" w x 28" d x 26" h 1,040mm x 712mm x 660mm h
*Outside Dimensions	55" w x 35" d x 38" h 1,400mm x 890mm x 965mm h
Approximate Volume	17.3 Cubic Feet 489 Liters
Electrical Requirements	N. American 115-120 Volts,60 Hz. 10 Amps Export Model 220-240 Volts,50 Hz. 5 Amps.
Maximum Working Pressure	6.0" of W.C. (1.49kPa)
Approximate Shipping (crated) Weight	250 Pounds 112.5 Kilos

*Includes Transfer Chamber

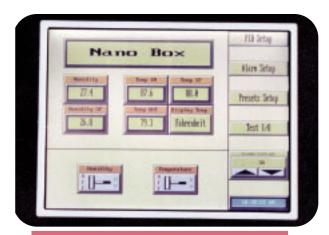

Close up of computer interface screen