

PURELAB® Dispenser

Uninterrupted Discovery



Proven reliability

Unrivalled purity design, unmatched flexibility & effortless use. Together with PURELAB Chorus & Quest, the PURELAB Dispenser enables uninterrupted discovery.



Accuracy:

Flexible dispensing with unrivalled built in purity monitoring sensor

- Delivering proof of water purity
- Repeated volume dispensing
- Improving lab efficiency
- Reproducible results

Streamlined footprint:

- - Maximized lab space
- - Optimization of storage space
- - Up to 4 dispensers per water purifier

Ultimate flexibility



Simple use and maintenance:

- Plug and Play installation with PURELAB range
- Height adjustable, 180 rotating arm & flexible dispenser arm
- IPX7 rating dispenser handset with
- Menu navigation
- Optional connection to point of use filtration

Efficient

- From precise drop by drop up to 2 litres per minute allowing for the right speed including calibration work
- Variable flow rate for easier filling of different sizes of labware
- Volumetric dispensing freeing you to continue working while water is being dispensed
- Profile & locked dispensing allowing you to choose your desired amount of water and prevent accidental flooding

*Dependent on water purifier purchased.
PURESURE is a twin-bed ion exchange process which allows the retention of any impurity released during exhaustion of the primary cartridge by the second polishing cartridge.
This process guarantees optimum product water quality, advanced warning of consumable change & extended consumable service life

*Maximum 5 meters away from the water purification system

Product Specifications

Height
580mm

Depth
375mm

Width
264mm

Intuitive dispenser handset

Multiple dispensing options

Point of Use filter

Height adjustable button

Flexible dispensing arm

Dispatch weight: 7.7kg

Operational Weight: 4.5kg

Installation: Free-standing on a worktop or secured to a worktop using the mounting bracket provided (see quick start guide for details).

Mains input

100-240 VAC, 50-60 Hz

Power required

24 V DC

Power consumption

120 VA

Noise Output

dBA - < 40

Inlet

8mm OD Tube


Return

8mm OD Tube

Flow bends

8mm

(provided in the basic installation kit with the main product)

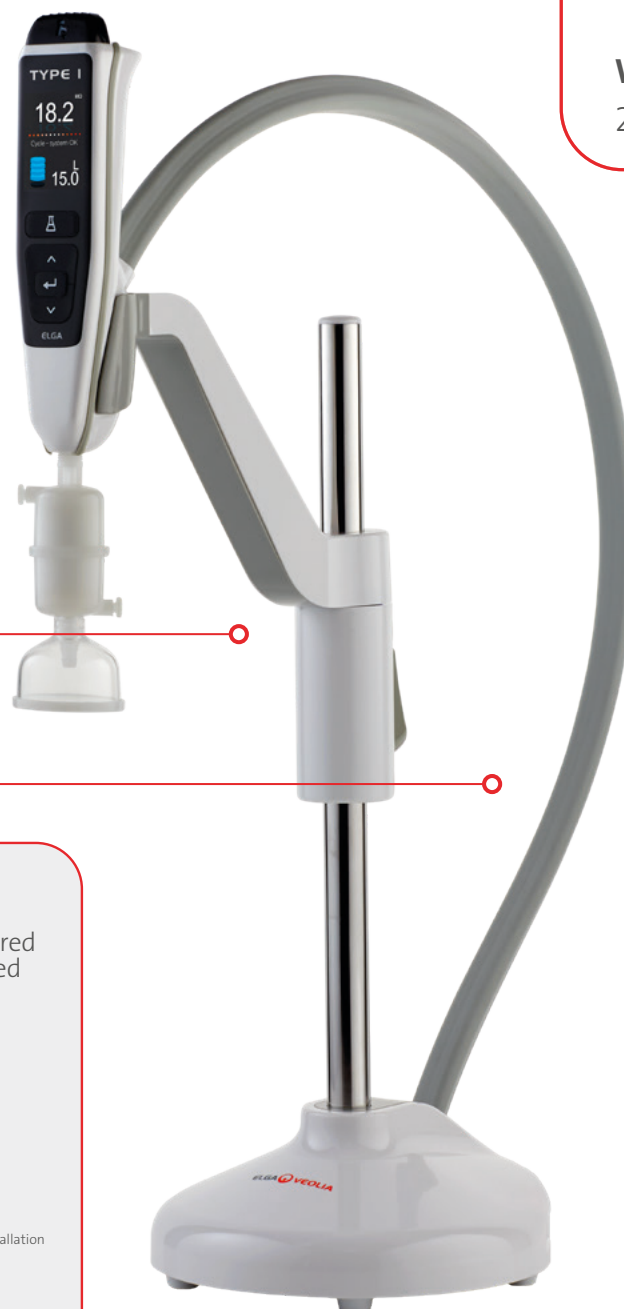

Dedicated to Discovery

ELGA Labwater are specialists in the engineering, service & support of water purification systems.

Unrivalled product design has achieved international recognition and awards.

Worldwide technical service teams support science & healthcare globally with specialist expertise.

Global digital performance monitoring from ensures laboratory work is uninterrupted.

A global supply chain supports clients from regional centres on every continent.

info@elgalabwater.com / www.elgalabwater.com

To find your nearest ELGA representative, go to www.elgalabwater.com and select your country for contact details.

Elga Global Operations Centre.

tel: +44 (0) 203 567 7300

fax: +44 (0) 203 567 7205